



DuPont™ AmberChrom™ XT20 Chromatography Resin


Reverse Phase Polymeric Resin for Purification and Polishing of Proteins, Peptides, and Oligonucleotides

Key Features

- Mechanical stability and chemical robustness to standard reversed-phase solvents and cleaning agents.
- Stable within a very broad pH range (up to pH 14), and easily cleaned in place (CIP) with most organic solvents and dilute acids and bases.
- Suitable for high-resolution, high-pressure chromatography.

Key Applications

- Capture, separation, and purification of peptides.
- Purification of oligonucleotides from impurities (failure sequences, DMT-off, oxidation products).
- Final polishing to remove trace impurities.

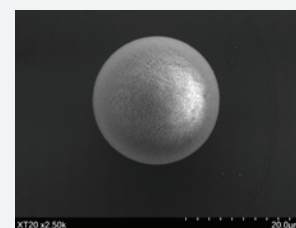

Typical Properties

Physical Properties

Copolymer	Crosslinked divinylbenzene
Matrix	Macroporous
Type	Adsorbent
Physical Form	White, opaque, spherical beads

Nitrogen BET

Surface Area	550 – 600 m ² /g
Total Pore Volume	0.61 mL/mL
Average Pore Diameter	300 Å

Chemical Properties

Functional Group	None
Shipping Form	Dry
Chemical Resistance	Insoluble in dilute solutions of acids or bases and common solvents: IPA, ACN, MeOH

Particle Size

Particle Diameter, mean	20 µm
15 – 30 µm	≥ 80%

Suggested Operating Conditions

Maximum Recommended Operating Temperature	60°C (140°F)
pH Range	1 – 14
Maximum Recommended Operating Pressure	60 bar (870 psi)

General Information

- DuPont™ AmberChrom™ XT20 chromatography resins can be used with medium pressure hardware (5 to 20 bar; 150 to 300 psi; 1000 kPa to 2000 kPa) or high pressure (HPLC column) up to 60 bar (880 psi/6000 kPa).
- DuPont™ AmberChrom™ XT20 chromatography resins are supplied in dry form and have a hydration ratio of 0.24 g/mL or 240 g/L. A slurry concentration of 40-45% is recommended for optimal column packing results, but higher percentage could be used if necessary.
- Store the column or used bulk resin in 20% solvent (propanol, methanol, ethanol, or acetonitrile), preferably between 4 – 25°C.

Important Information

- Polymeric adsorbents, as produced, contain by-products resulting from the manufacturing process. The user must determine the extent to which organic by-product must be removed for any particular use and establish techniques to assure that the appropriate level of purity is achieved for that use.
- Like any chromatographic stationary phase, a conditioning step with the working solvent must be performed before operation.
- **WARNING:** Oxidizing agents such as nitric acid attack organic ion exchange resins under certain conditions. This could lead to anything from slight resin degradation to a violent exothermic reaction (explosion). Before using strong oxidizing agents, consult sources knowledgeable in handling such materials.

Regulatory Note

DuPont can provide regulatory support for DuPont™ AmberChrom™ XT20 chromatography resins to end users under confidentiality, upon request.



Have a question? Contact us at:
[dupont.com/water/contact-us](https://www.dupont.com/water/contact-us)

All information set forth herein is for informational purposes only. This information is general information and may differ from that based on actual conditions. Customer is responsible for determining whether products and the information in this document are appropriate for Customer's use and for ensuring that Customer's workplace and disposal practices are in compliance with applicable laws and other government enactments. The product shown in this literature may not be available for sale and/or available in all geographies where DuPont is represented. The claims made may not have been approved for use in all countries. Please note that physical properties may vary depending on certain conditions and while operating conditions stated in this document are intended to lengthen product lifespan and/or improve product performance, it will ultimately depend on actual circumstances and is in no event a guarantee of achieving any specific results. DUPONT ASSUMES NO OBLIGATION OR LIABILITY FOR THE INFORMATION IN THIS DOCUMENT. References to "DuPont" or the "Company" mean the DuPont legal entity selling the products to Customer unless otherwise expressly noted. NO WARRANTIES ARE GIVEN; ALL IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE ARE EXPRESSLY EXCLUDED. No freedom from infringement of any patent or trademark owned by DuPont or others is to be inferred.

DuPont™, the DuPont Oval Logo, and all trademarks and service marks denoted with ™, SM or ® are owned by affiliates of DuPont de Nemours Inc. unless otherwise noted. © 2025 DuPont. All rights reserved.

Form No. 45-D04054-en, Rev. 1
April 2025